


Ultimate Formula to Monitor Diabetes

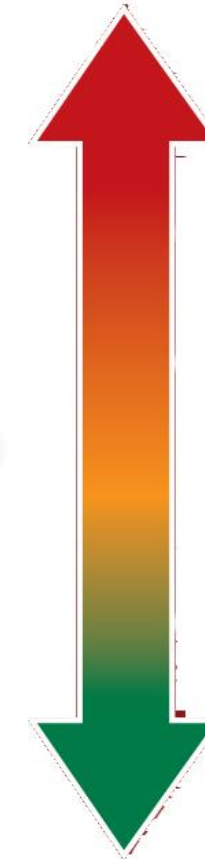

Action Suggested	HbA1c Test Score	Estimated Average Glucose mg/dl	mmol/L
	Excellent	4.0	50
Good	5.0	80	4.7
	6.0	115	6.3
	7.0	150	8.2
Action Suggested	8.0	180	10.0
	9.0	215	11.9
	10.0	250	13.7
	11.0	280	15.6
	12.0	315	17.4
	13.0	350	19.3
	14.0	380	21.1

HbA1c and eAG in one sample

Intelligent calibration using radio frequency technology

- No interference from Hb Variants
- Easy access to PC
- Accurate result in 7 minutes
- Cost effectiveness per test
- Patient archiving software
- Portable with built-in battery


Twin A1C Test



Blood Sample Preparation



Hemolysate Preparation

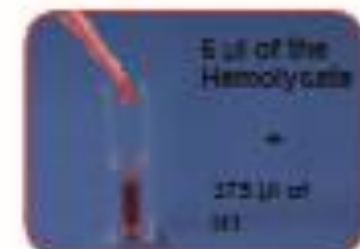

Mix Carefully



Locate in reading location & press enter



After 2 min add 75 µL of R2



Record the result